



Semester: I

CHEMISTRY OF SMART MATERIALS AND DEVICES (22CHY12A)

(Category: Professional Core Course) Stream: CS (Theory and Practice)

Practical:

Expt-1: Estimation of Copper from PCBs e-waste

C. MANJUNATHA, Ph.D, MRSC

Assistant Professor, Department of Chemistry, RV College of Engineering

Bengaluru-560059, INDIA, Phone: +91 9036651277

E-mail: manjunathac@rvce.edu.in; manju.chem20@gmail.com

Estimation of Copper from PCBs e-waste



“E-waste”, “electronic waste”, “e-scrap” Are used electronics that are nearing the end of their useful life, and are discarded.


Waste printed circuit boards (PCBs)

- ❑ PCBs comprise metals (Cu, Sn, Pb, Ag, Au, Pd, Fe, Ni, and Cr), nonmetals (glass fibers, electronic component insulators, capacitors, resistors, and so on), and organic compounds (epoxy resin, paints...).
- ❑ Copper is the primary component of PCB, which is employed as an electric current conductor.
- ❑ These discarded metals without treatment threaten the economy, the environment, and human health.

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)




1. All the detachable components free from copper (e.g., relays, capacitors, and resistance) were separated, the bare PCB boards were shredded and crushed by a cutting and grinding machine and passed through a sieve mesh to obtain the WPCB powder of varying particle sizes.
2. Then the mixer of WPCB powder and $K_2S_2O_8$ (with suitable mass ratios) together with grinding balls, were sealed in the ZrO_2 pot. The mixture was then co-ground at different rotary speeds under ambient conditions for different periods of time (**Stage I**).
3. Co-ground products of PCB powder and $K_2S_2O_8$ is leached in deionized water with constant magnetic stirring at different temperatures: 25, 35, 45, and 55 °C. The leaching solution and solid residues were rapidly separated by vacuum filtration after the reaction was complete. (**Stage-II**).
4. After filtration, $CuSO_4 \cdot 5H_2O$ could be easily separated from K_2SO_4 solution because of their differences in solubility. (**Stage-III**)

The flowchart for Cu recovery from PCB powder is shown in Figure.

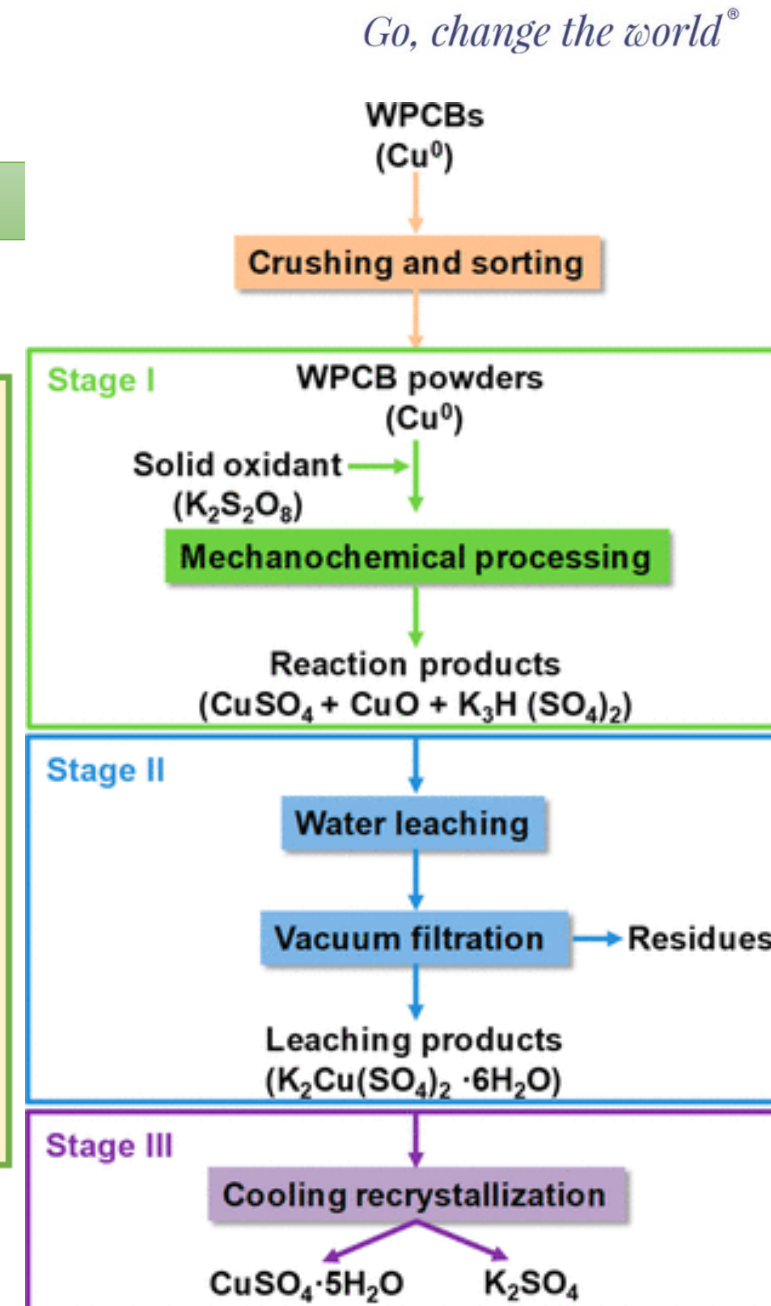

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

How Cu is
estimated?
Iodometric
method.

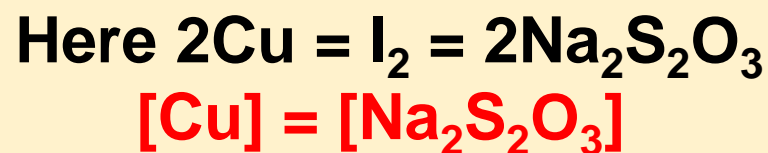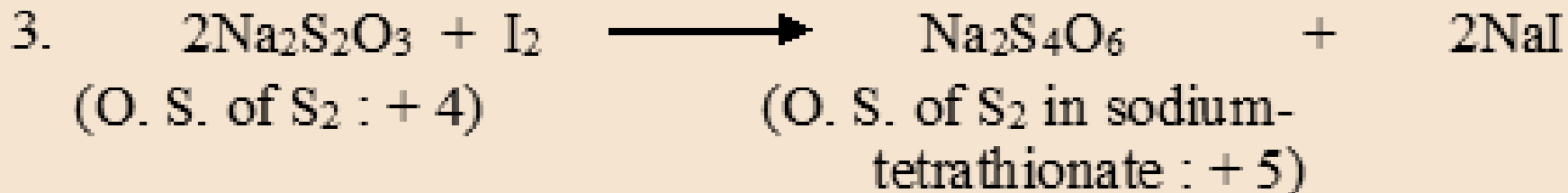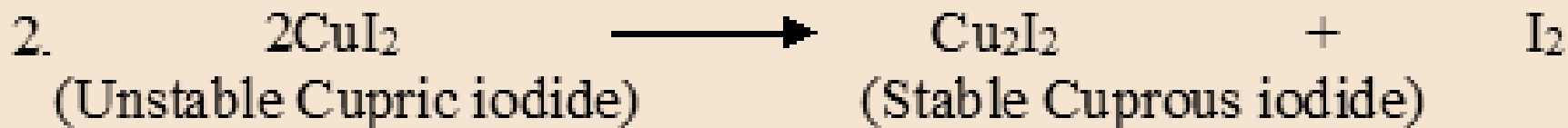

❖ The I_2 liberated is estimated using Std. $\text{Na}_2\text{S}_2\text{O}_3$ solution by **IODOMETRIC** titration method.

❖ **The amount of $\text{Na}_2\text{S}_2\text{O}_3$ consumed is equivalent to $[\text{Cu}]$ present in water.**

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

Estimation of Copper from PCBs e-waste

Model Procedure:

Step-1 : Extraction of copper from e-waste (PCBs)

After removing all the detachable components, the PCB boards are crushed to powder is mixed with $K_2S_2O_8$ (Oxidant) and ball milled to fine powder. The fine mixer of PCB+ $K_2S_2O_8$ powder is leached in deionized water and nitric acid with constant stirring at 55 °C. The leached solution and solid residue is separated by filtration. After the filtration, the copper solution is mixed with nitric acid to prevent hydrolysis. (called PCB solution)

Step-2: Preparation of Standard Sodium Thiosulfate solution ($Na_2S_2O_3$)

The accurate quantity ('w' g) of $Na_2S_2O_3 \cdot 5H_2O$ (Molecular Weight =Equivalent wt=248.18) is weighed and dissolved in 100 ml of pure distilled water in 100 ml volumetric flask.

$$\text{Normality of } Na_2S_2O_3 \cdot 5H_2O = (w \times 10) / \text{Equivalent wt} = 'Y' N$$

Step-3 Titration

Burette: Standard ('X'N) $Na_2S_2O_3$ solution

Conical flask: 25 ml of PCB solution + 1 spatula urea + Heat near to boil + cool to room temperature + Ammonium hydroxide drop by drop till to get bluish with precipitate + $\frac{1}{4}$ test tube acetic Acid (~5 ml) + 1 test tube of KI (~15 ml)(dark brownish yellow solution with white precipitate appears) + titrate against $Na_2S_2O_3$ solution till to get pale yellow color + add 5 drops of starch indicator+ solution turns blue+ continue titration till to get permanent pure white solution. (End point Reached)

Note down the volume of 'V' ml for estimation of copper in PCB solution.


$$1000 \text{ ml of PCB solution contains} = \{V \times Y \times 63.54\} / 25 \text{ g}$$

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

Role of chemicals added



Urea :

It removes nitrogen oxides

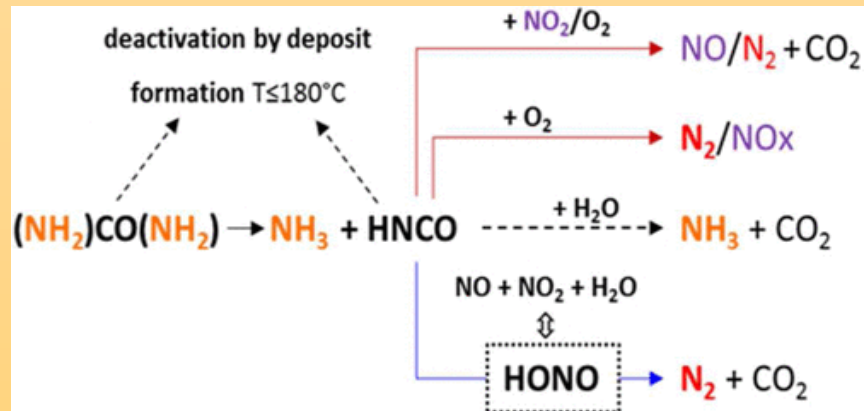

Ammonium hydroxide :

It neutralises the excess H^+ ions (Acidic to basic)

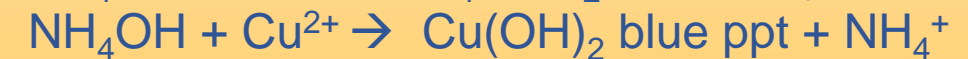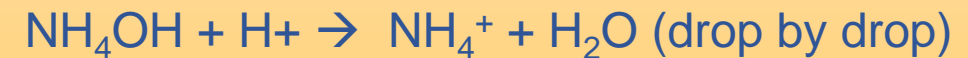

Acetic acid:

It neutralises the $\text{Cu}(\text{OH})_2$ and provides the required slightly acidic medium of $\text{pH} = 3-4$.





Potassium Iodide (KI):

It is a special chemical, which reduces only copper ions from Cu^{2+} to Cu^+ , in the presence of other metal ions in PCB solution. And liberate equivalent amount of Iodine (I_2)

(Reaction given in previous slide)

Std. $\text{Na}_2\text{S}_2\text{O}_3$
(‘Y’N)

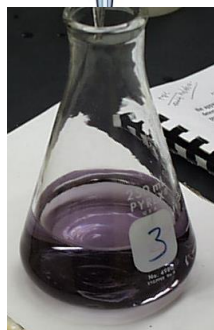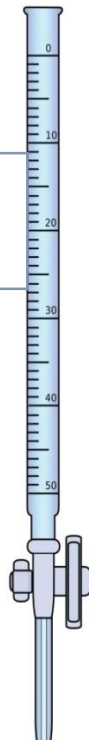


Procedure : Iodometric method.

The end point
of the titration
is BLUE to
COLORLESS


25ml of PCB solution+ One spatula
Urea+ Ammonium hydroxide d/d +
1/2 test tube of Acetic acid + 1 test
tube of KI + Starch indicator.

Step (1): Reacting $[\text{Cu}^{2+}]$ with KI and liberating I_2

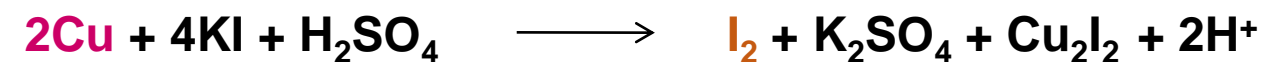

The amount of I_2 liberated is equivalent to $[\text{Cu}]$.

Step (2): Reacting I_2 with Std $\text{Na}_2\text{S}_2\text{O}_3$

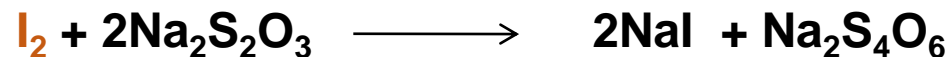


Let the end point of the titration is = ‘V’ ml. volume of the

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

When Iodine (DARK Brown) is liberated, which is equivalent to Cu

When it is titrated with $\text{Na}_2\text{S}_2\text{O}_3$, Iodine reacts and its concentration decreases (Brown color decreases)

When starch indicator is added, Iodine reacts with starch, forms dark blue (violet) solution.

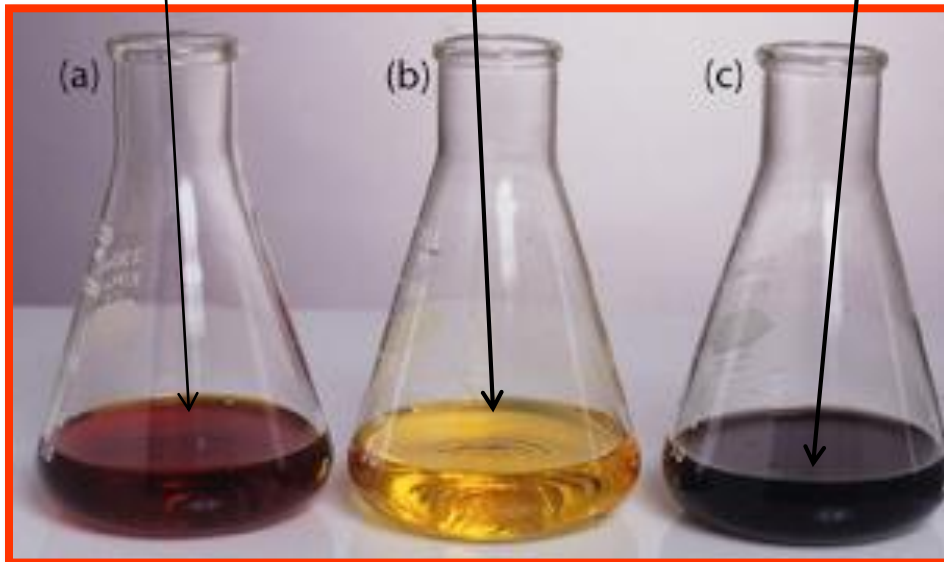


Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

AMYLOSE

Stains deep blue with iodine
Relative molecular mass up to 50 000
Up to 300 glucose units/molecules
Unbranched helical chain

STRUCTURE OF MOLECULE

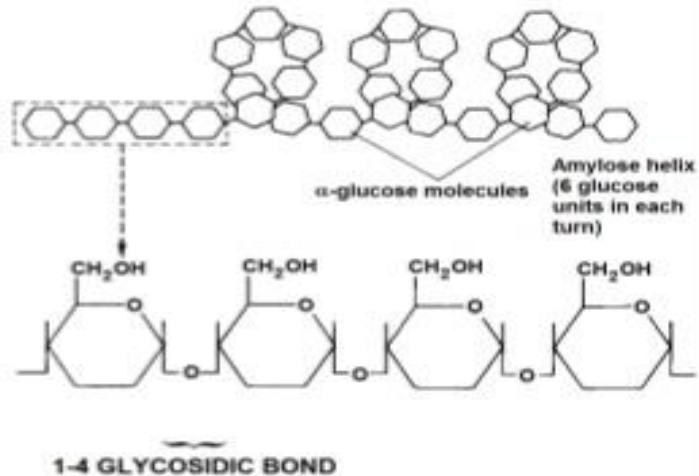


**Linear
molecule**

AMYLOPECTIN

Stains red to purple with iodine
Relative molecular mass up to 500 000
1300-1500 glucose units/molecules
Branched chain

STRUCTURE OF MOLECULE

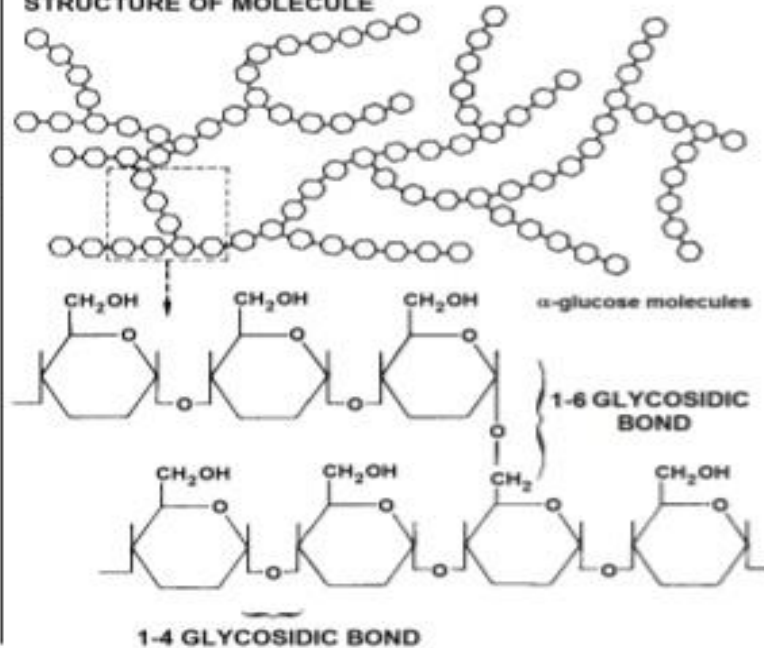


Branched molecule



Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

How it works – Mechanism of Reaction

3. The free iodine combines with iodide very rapidly to form the negative ion I_3^- , which reacts with starch to form a dark blue complex:

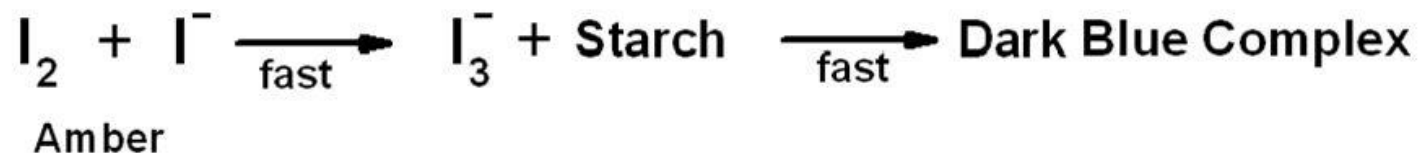

Starch - Iodine Complex

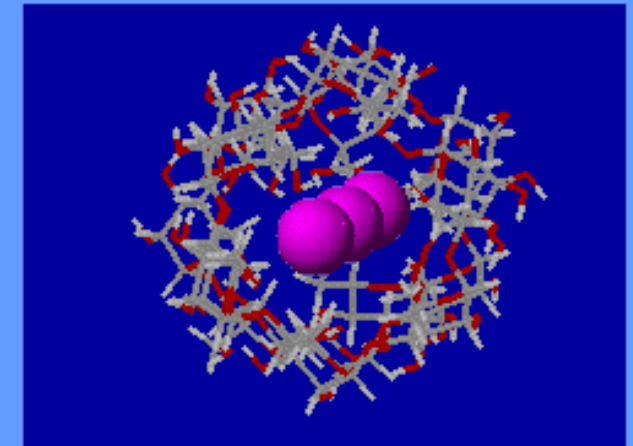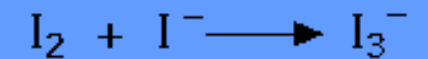

Iodine slides into starch coil to give a blue-black color

Estimation of Copper from PCBs e-waste (Waste PCBs: WPCB)

Observations and calculations:

Note: The equivalent wt of Cu is 63.54 g (The change in oxidation number is 1)

One gram equivalent wt of $\text{Na}_2\text{S}_2\text{O}_3$ = one gram equivalent wt of Copper

i.e, 1000 ml of 1N $\text{Na}_2\text{S}_2\text{O}_3$ = 63.54 g (Eq. Wt.) of Copper

Therefore, ' V ' ml of ' Y ' N $\text{Na}_2\text{S}_2\text{O}_3$ = ?????

= ($V \times Y \times 63.54$) / 1000 in grams

i.e, 25 ml of PCB solution titrated contains = ($V \times Y \times 63.54$) / 1000 in grams

Therefore, 1000 ml of PCB solution contains= ????

= { $V \times Y \times 63.54$ } / 25

=g/L

The amount of Cu the given 1000 ml of PCB solution is =g